

Load Balancing Pathomation PMA.core & PMA.view

Version 1.0.0



Table of Contents

1. About this Brief	3
2. Loadbalancer.org Appliances Supported	3
3. Software Versions Supported	3
3.1. Loadbalancer.org Appliance	3
3.2. Pathomation PMA.core & PMA.view	3
4. Pathomation Platform	3
5. Load Balancing PMA.core & PMA.view	3
5.1. Load Balancing & HA Requirements	3
5.2. Virtual Service (VIP) Requirements	4
6. Deployment Concept	4
7. Load Balancer Deployment Methods	4
7.1. Layer 7 SNAT Mode	5
8. Loadbalancer.org Appliance – the Basics	6
8.1. Virtual Appliance	
8.2. Initial Network Configuration	6
8.3. Accessing the Appliance WebUI	6
8.3.1. Main Menu Options	8
8.4. Appliance Software Update.	
8.4.1. Online Update	
8.4.2. Offline Update	
8.5. Ports Used by the Appliance	
8.6. HA Clustered Pair Configuration	
9. Appliance Configuration for Pathomation PMA.core & PMA.view	
9.1. VIP 1 - PMA.core	
9.1.1. Virtual Service (VIP) Configuration.	
9.1.2. Configure the Associated Real Servers (RIPs)	
9.2. VIP 2 - PMA.view	
9.2.1. Virtual Service (VIP) Configuration.	
9.2.2. Configure the Associated Real Servers (RIPs)	
9.3. Finalizing the Configuration	
10. Testing & Verification	
10.1. Accessing PMA.core & PMA.view via the Load Balancer	
10.2. Using System Overview	
11. Technical Support	
12. Further Documentation	
13. Appendix	
13.1. Configuring HA - Adding a Secondary Appliance	
13.1.1. Non-Replicated Settings	
13.1.2. Configuring the HA Clustered Pair	
14. Document Revision History	18

1. About this Brief

This brief outlines the steps required to configure a load balanced Pathomation PMA.core & PMA.view environment utilizing Loadbalancer.org appliances. It covers the configuration of the load balancers and also any Pathomation platform configuration changes that are required to enable load balancing.

For more information about initial appliance deployment, network configuration and using the Web User Interface (WebUI), please also refer to the Administration Manual.

2. Loadbalancer.org Appliances Supported

All our products can be used with PMA.core & PMA.view. For full specifications of available models, please refer to https://www.loadbalancer.org/products/enterprise.

Some features may not be available or fully supported in all cloud platforms due to platform specific limitations. For more details, please refer to the "Main Differences to our Standard (Non-Cloud) Product" section in the appropriate cloud platform Quick Start Guide or check with Loadbalancer.org support.

3. Software Versions Supported

3.1. Loadbalancer.org Appliance

v8.13.2 and later

8 Note

The screenshots used throughout this document aim to track the latest Loadbalancer.org software version. If you're using an older version, or the very latest, the screenshots presented here may not match your WebUI exactly.

3.2. Pathomation PMA.core & PMA.view

• v3.0 and later

4. Pathomation Platform

PMA.core and PMA.view are components of the Pathomation platform. The platform is used for managing, distributing, viewing and annotating digital images of human tissue slides, scanned using any commercially available slide scanner capable of producing whole slide images in a supported file format.

5. Load Balancing PMA.core & PMA.view

8 Note

It's highly recommended that you have a working Pathomation Platform environment first before implementing the load balancer.

5.1. Load Balancing & HA Requirements

Both PMA.core and PMA.view can be load balanced across multiple servers to ensure high availability. Typically,



PMA.core and PMA.view are installed on each physical server and accessed on the same ports (80,443), but through different URL paths (e.g., http://server_ip/PMA.core and http://server_ip/PMA.view).

Separate Virtual Servers (VIPs) can be configured for each component (PMA.core and PMA.view) which allows independent load balancing and health checks for each application.

Shared backend resources, including the SQL Server database and image storage, ensure that all nodes provide consistent data access and session persistence for web-based image viewing and annotation. Proper configuration of Host headers or URL paths may be required to direct traffic correctly to each component.


5.2. Virtual Service (VIP) Requirements

To provide load balancing and HA for Pathomation PMA.core & PMA.view, the following VIPs are required:

Ref.	VIP Name	Mode	Port(s)	Persistence Mode	Health Check
VIP 1	PMA.core	L7 SNAT	80,443	Source IP	Negotiate HTTP (GET)
VIP 2	PMA.view	L7 SNAT	80,443	Source IP	Negotiate HTTP (GET)

6. Deployment Concept

Once the load balancer is deployed, clients connect to the Virtual Service (VIP) rather than connecting directly to one of the Pathomation PMA.core and PMA.view servers. These connections are then load balanced across the PMA.core and PMA.view servers to distribute the load according to the load balancing algorithm selected.

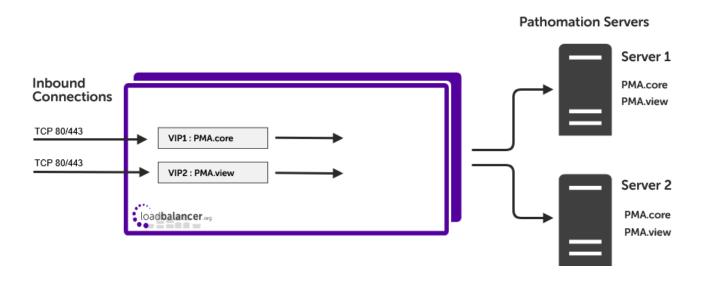

8 Note

The load balancer can be deployed as a single unit, although Loadbalancer.org recommends a clustered pair for resilience & high availability. Please refer to the section Configuring HA - Adding a Secondary Appliance in the appendix for more details on configuring a clustered pair.

7. Load Balancer Deployment Methods

The load balancer can be deployed in 4 fundamental ways: Layer 4 DR mode, Layer 4 NAT mode, Layer 4 SNAT




mode, and Layer 7 SNAT mode.

For Pathomation PMA.core & PMA.view, layer 7 SNAT mode is recommended. This mode is described below and is used for the configuration presented in this guide.

7.1. Layer 7 SNAT Mode

Layer 7 SNAT mode uses a proxy (HAProxy) at the application layer. Inbound requests are terminated on the load balancer and HAProxy generates a new corresponding request to the chosen Real Server. As a result, Layer 7 is typically not as fast as the Layer 4 methods. Layer 7 is typically chosen when either enhanced options such as SSL termination, cookie based persistence, URL rewriting, header insertion/deletion etc. are required, or when the network topology prohibits the use of the layer 4 methods. The image below shows an example network diagram for this mode.

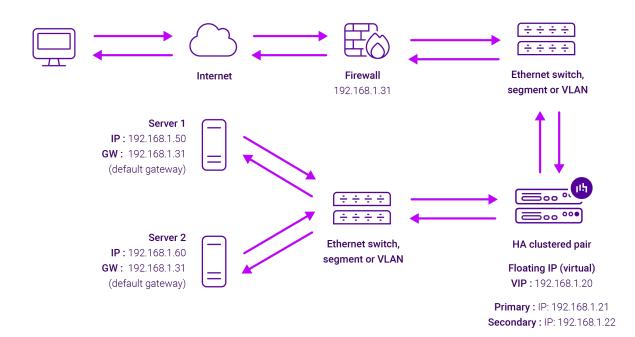

- Because layer 7 SNAT mode is a full proxy, Real Servers in the cluster can be on any accessible network including across the Internet or WAN.
- Layer 7 SNAT mode is not transparent by default, i.e. the Real Servers will not see the source IP address of the client, they will see the load balancer's own IP address by default, or any other local appliance IP address if preferred (e.g. the VIP address). This can be configured per layer 7 VIP. If required, the load balancer can be configured to provide the actual client IP address to the Real Servers in 2 ways. Either by inserting a header that contains the client's source IP address, or by modifying the Source Address field of the IP packets and replacing the IP address of the load balancer with the IP address of the client. For more information on these methods please refer to Transparency at Layer 7.
- Layer 7 SNAT mode can be deployed using either a one-arm or two-arm configuration. For two-arm
 deployments, eth1 is typically used for client side connections and eth0 is used for Real Server connections,
 although this is not mandatory since any interface can be used for any purpose.
- Requires no mode-specific configuration changes to the load balanced Real Servers.
- Port translation is possible with Layer 7 SNAT mode, e.g. VIP:80 → RIP:8080 is supported.

• You should not use the same RIP:PORT combination for layer 7 SNAT mode VIPs and layer 4 SNAT mode VIPs because the required firewall rules conflict.

8. Loadbalancer.org Appliance – the Basics

8.1. Virtual Appliance

A fully featured, fully supported 30 day trial is available if you are conducting a PoC (Proof of Concept) deployment. The VA is currently available for VMware, Virtual Box, Hyper-V, KVM, XEN and Nutanix AHV and has been optimized for each Hypervisor. By default, the VA is allocated 2 vCPUs, 4GB of RAM and has a 20GB virtual disk. The Virtual Appliance can be downloaded here.

8 Note	The same download is used for the licensed product, the only difference is that a license key file (supplied by our sales team when the product is purchased) must be applied using the appliance's WebUI.
8 Note	Please refer to Virtual Appliance Installation and the ReadMe.txt text file included in the VA download for additional information on deploying the VA using the various Hypervisors.
	download for additional information on deploying the vA using the various hypervisors.
ੀ Note	The VA has 4 network adapters. For VMware only the first adapter (eth0) is connected by default. For HyperV, KVM, XEN and Nutanix AHV all adapters are disconnected by default. Use the network configuration screen within the Hypervisor to connect the required adapters.

8.2. Initial Network Configuration

After boot up, follow the instructions on the appliance console to configure the management IP address, subnet mask, default gateway, DNS servers and other network and administrative settings.

(!) Important Be sure to set a secure password for the load balancer, when prompted during the setup routine.

8.3. Accessing the Appliance WebUI

The WebUI is accessed using a web browser. By default, users are authenticated using Apache authentication. Users can also be authenticated against LDAP, LDAPS, Active Directory or Radius - for more information, please refer to External Authentication.

Note

There are certain differences when accessing the WebUI for the cloud appliances. For details, please refer to the relevant Quick Start / Configuration Guide.

1. Using a browser, navigate to the following URL:

https://<IP-address-configured-during-the-network-setup-wizard>:9443/lbadmin/

Note
You'll receive a warning about the WebUl's SSL certificate. This is due to the default self signed certificate that is used. If preferred, you can upload your own certificate - for more

information, please refer to Appliance Security Features.

8 Note

If you need to change the port, IP address or protocol that the WebUI listens on, please refer to Service Socket Addresses.

2. Log in to the WebUI using the following credentials:

Username: loadbalancer

Password: <configured-during-network-setup-wizard>


8 Note

To change the password, use the WebUI menu option: *Maintenance > Passwords*.

Once logged in, the WebUI will be displayed as shown below:





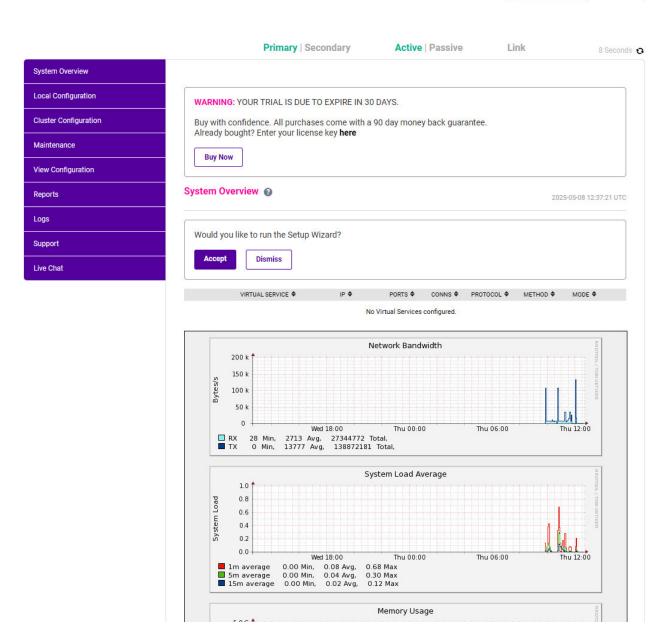

3. You'll be asked if you want to run the Setup Wizard. Click **Dismiss** if you're following a guide or want to configure the appliance manually. Click **Accept** to start the Setup Wizard.

8 Note

The Setup Wizard can only be used to configure Layer 7 services.

8.3.1. Main Menu Options

System Overview - Displays a graphical summary of all VIPs, RIPs and key appliance statistics

Local Configuration - Configure local host settings such as IP address, DNS, system time etc.

Cluster Configuration - Configure load balanced services such as VIPs & RIPs

Maintenance - Perform maintenance tasks such as service restarts and creating backups

View Configuration - Display the saved appliance configuration settings

Reports - View various appliance reports & graphs

Logs - View various appliance logs

Support - Create a support download, contact the support team & access useful links

Live Chat - Start a live chat session with one of our Support Engineers

8.4. Appliance Software Update

We recommend that the appliance is kept up to date to ensure that you benefit from the latest bug fixes, security updates and feature improvements. Both online and offline update are supported.

Note For full details, please refer to Appliance Software Update in the Administration Manual.

8 Note

Services may need to be restarted/reloaded after the update process completes or in some cases a full appliance restart may be required. We therefore recommend performing the update during a maintenance window.

8.4.1. Online Update

The appliance periodically contacts the Loadbalancer.org update server (**update.loadbalancer.org**) and checks for updates. This is the default behavior and can be disabled if preferred. If an update is found, a notification similar to the example below will be displayed at the top of the WebUI:

Information: Update 8.13.2 is now available for this appliance.

Online Update

Click **Online Update**. A summary of all new features, improvements, bug fixes and security updates included in the update will be displayed. Click **Update** at the bottom of the page to start the update process.

(1) Important Do not navigate away whilst the update is ongoing, this may cause the update to fail.

The update can take several minutes depending on download speed and upgrade version. Once complete, the following message will be displayed:



Information: Update completed successfully. Return to system overview.

If services need to be reloaded/restarted or the appliance needs a full restart, you'll be prompted accordingly.

8.4.2. Offline Update

If the appliance does not have access to the Internet, offline update can be used.

To check for the latest version, please refer to our product roadmap page available here. To obtain the latest offline update files contact support@loadbalancer.org.

To perform an offline update:


- 1. Using the WebUI, navigate to: *Maintenance > Software Update*.
- 2. Select Offline Update.
- 3. The following screen will be displayed:

Software Update

Offline Update

The following steps will lead you through offline update.

- 1. Contact Loadbalancer.org support to obtain the offline update archive and checksum.
- 2. Save the archive and checksum to your local machine.
- 3. Select the archive and checksum files in the upload form below.
- 4. Click Upload and Install to begin the update process.

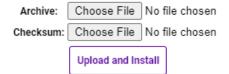

- 4. Select the Archive and Checksum files.
- 5. Click Upload and Install.
- 6. If services need to be reloaded/restarted or the appliance needs a full restart, you'll be prompted accordingly.

8.5. Ports Used by the Appliance

By default, the appliance uses the following TCP & UDP ports:

Protocol	Port	Purpose
TCP	22 *	SSH
TCP & UDP	53 *	DNS / GSLB
TCP & UDP	123	NTP



Protocol	Port	Purpose
TCP & UDP	161 *	SNMP
UDP	6694	Heartbeat between Primary & Secondary appliances in HA mode
TCP	7778	HAProxy persistence table replication
TCP	9000 *	Gateway service (Centralized/Portal Management)
TCP	9080 *	WebUI - HTTP (disabled by default)
TCP	9081 *	Nginx fallback page
TCP	9443 *	WebUI - HTTPS
TCP	25565 *	Shuttle service (Centralized/Portal Management)

The ports used for SSH, GSLB, SNMP, the WebUI, the fallback page, the gateway service and the shuttle service can be changed if required. For more information, please refer to Service Socket Addresses.

8.6. HA Clustered Pair Configuration


Loadbalancer.org recommend that load balancer appliances are deployed in pairs for high availability. In this guide a single unit is deployed first, adding a secondary unit is covered in the section Configuring HA - Adding a Secondary Appliance of the appendix.

9. Appliance Configuration for Pathomation PMA.core & PMA.view

9.1. VIP 1 - PMA.core

9.1.1. Virtual Service (VIP) Configuration


1. Using the WebUI, navigate to *Cluster Configuration > Layer 7 - Virtual Services* and click **Add a new Virtual Service**.

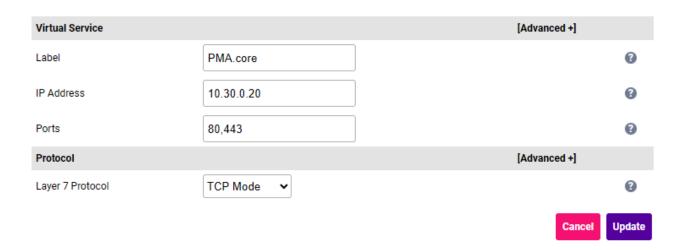

- 2. Enter a suitable Label (name) for the Virtual Service, e.g. PMA.core.
- 3. Set the Virtual Service IP Address field to the required IP address, e.g. 10.30.0.20.
- 4. Set Ports to 80,443.
- 5. Set the Layer 7 Protocol to TCP Mode.
- 6. Click Update.
- 7. Now click **Modify** next to the newly created VIP.
- 8. Scroll down to the *Health Checks* section.
 - Set the Health Checks to Negotiate HTTP (GET).
- 9. Click Update.

9.1.2. Configure the Associated Real Servers (RIPs)

1. Using the WebUI, navigate to Cluster Configuration > Layer 7 - Real Servers and click on Add a new Real Server next to the newly created VIP.

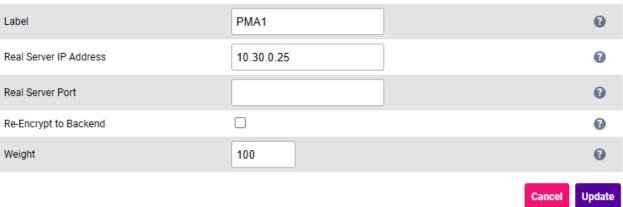




- 2. Enter a suitable Label (name) for the Real Server, e.g. PMA1.
- 3. Set the Real Server IP Address field to the required IP address, e.g. 10.30.0.25.
- 4. Leave the Real Server Port field blank.
- 5. Click Update.
- 6. Repeat these steps to add additional Real Servers as required.

9.2. VIP 2 - PMA.view

9.2.1. Virtual Service (VIP) Configuration


1. Using the WebUI, navigate to Cluster Configuration > Layer 7 - Virtual Services and click Add a new Virtual Service.

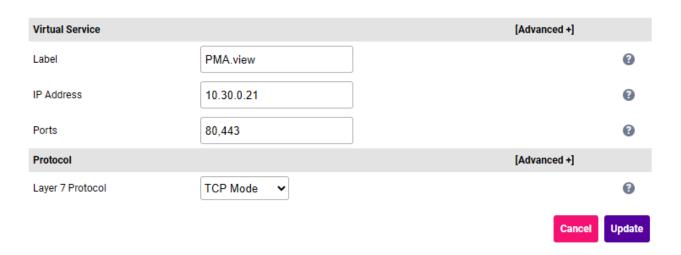

- 2. Enter a suitable Label (name) for the Virtual Service, e.g. PMA.view.
- 3. Set the Virtual Service IP Address field to the required IP address, e.g. 10.30.0.21.
- 4. Set Ports to 80,443.
- 5. Set the *Layer 7 Protocol* to **TCP Mode**.
- 6. Click Update.
- 7. Now click **Modify** next to the newly created VIP.
- 8. Scroll down to the *Health Checks* section.
 - Set the *Health Checks* to **Negotiate HTTP (GET)**.
- 9. Click Update.

9.2.2. Configure the Associated Real Servers (RIPs)

1. Using the WebUI, navigate to *Cluster Configuration > Layer 7 – Real Servers* and click on **Add a new Real Server** next to the newly created VIP.

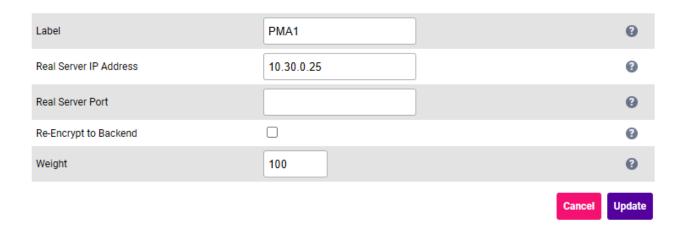

- 2. Enter a suitable Label (name) for the Real Server, e.g. PMA1.
- 3. Set the Real Server IP Address field to the required IP address, e.g. 10.30.0.25.
- 4. Leave the Real Server Port field blank.

- 5. Click Update.
- 6. Repeat these steps to add additional Real Servers as required.

9.3. Finalizing the Configuration

To apply the new settings, HAProxy must be reloaded. This can be done using the button in the "Commit changes" box at the top of the screen or by using the Restart Services menu option:

- 1. Using the WebUI, navigate to: Maintenance > Restart Services.
- 2. Click Reload HAProxy.

10. Testing & Verification

8 Note

For additional guidance on diagnosing and resolving any issues you may have, please also refer to Diagnostics & Troubleshooting.

10.1. Accessing PMA.core & PMA.view via the Load Balancer

Verify that you're able to successfully access PMA.core & PMA.view via the Virtual Services on the load balancer.

PMA.core

- 1. Basic Web Access Test
 - Open a browser and navigate to the PMA.core VIP (HTTP/HTTPS).
 - Confirm the login page loads without error.
- 2. Image Viewing Test
 - Open a whole-slide image from the web interface.
 - Pan/zoom around the slide and verify tiles load quickly without errors (tests tile streaming across multiple requests).
- 3. Upload Test
 - Upload a sample slide via the web interface.
 - Confirm it appears in the shared repository and is viewable from another client session.

PMA.view

- 1. Basic Web Access Test
 - Open a browser and navigate to the PMA.view VIP (HTTP/HTTPS).
 - Confirm the login page loads without error.
- 2. Image Viewing Test
 - Open a whole-slide image from the web interface.

 Pan/zoom around the slide and verify tiles load quickly without errors (tests tile streaming across multiple requests).

3. Slide Load Test

- Open a whole-slide image through PMA.view.
- Verify that tiles stream seamlessly with no missing or corrupted sections.

4. Pan & Zoom Test

- Pan, zoom in, and zoom out across the image.
- Confirm responsiveness and smoothness, ensuring the load balancer distributes traffic without latency spikes.

8 Note

Make sure that DNS is updated so that any FQDNs used point to the VIPs rather than individual servers.

10.2. Using System Overview

The System Overview can be viewed in the WebUI. It shows a graphical view of all Virtual Services & the associated Real Servers (i.e. the PMA.core & PMA.view servers) and shows the state/health of each server as well as the overall state of each cluster. The example below shows that all servers are healthy (green) and available to accept connections:

2025-10-21 10:53:33 UTC

VIRTUAL SERVICE ♦	IP ♦	PORTS ♦	CONNS ♦	PROTOCOL ♦	METHOD	♦ MODE ♦	
PMA.core	10.30.0.20	80,443	0	TCP	Layer 7	Proxy	9,4
REAL SERVER	IP	PORTS	WEIGHT	CONNS			
PMA1	10.30.0.25	80,443	100	0	Drain	Halt	W
PMA2	10.30.0.26	80,443	100	0	Drain	Halt	9,4
PMA.view	10.30.0.21	80,443	0	TCP	Layer 7	Proxy	W
REAL SERVER	IP	PORTS	WEIGHT	CONNS			
PMA1	10.30.0.25	80,443	100	0	Drain	Halt	9,1
							9.4
	PMA.core REAL SERVER PMA1 PMA2 PMA.view REAL SERVER	PMA.core 10.30.0.20 REAL SERVER IP PMA1 10.30.0.25 PMA2 10.30.0.26 PMA.view 10.30.0.21 REAL SERVER IP	PMA.core 10.30.0.20 80,443 REAL SERVER IP PORTS PMA1 10.30.0.25 80,443 PMA2 10.30.0.26 80,443 PMA.view 10.30.0.21 80,443 REAL SERVER IP PORTS	PMA.core 10.30.0.20 80,443 0 REAL SERVER IP PORTS WEIGHT PMA1 10.30.0.25 80,443 100 PMA2 10.30.0.26 80,443 100 PMA.view 10.30.0.21 80,443 0 REAL SERVER IP PORTS WEIGHT	PMA.core 10.30.0.20 80,443 0 TCP REAL SERVER IP PORTS WEIGHT CONNS PMA1 10.30.0.25 80,443 100 0 PMA2 10.30.0.26 80,443 100 0 PMA.view 10.30.0.21 80,443 0 TCP REAL SERVER IP PORTS WEIGHT CONNS	PMA.core 10.30.0.20 80,443 0 TCP Layer 7 REAL SERVER IP PORTS WEIGHT CONNS PMA1 10.30.0.25 80,443 100 0 Drain PMA2 10.30.0.26 80,443 100 0 Drain PMA.view 10.30.0.21 80,443 0 TCP Layer 7 REAL SERVER IP PORTS WEIGHT CONNS	PMA.core 10.30.0.20 80,443 0 TCP Layer 7 Proxy REAL SERVER IP PORTS WEIGHT CONNS PMA1 10.30.0.25 80,443 100 0 Drain Halt PMA2 10.30.0.26 80,443 100 0 Drain Halt PMA.view 10.30.0.21 80,443 0 TCP Layer 7 Proxy REAL SERVER IP PORTS WEIGHT CONNS

11. Technical Support

For more details about configuring the appliance and assistance with designing your deployment please don't hesitate to contact the support team using the following email address: support@loadbalancer.org.

12. Further Documentation

For additional information, please refer to the Administration Manual.



13. Appendix

13.1. Configuring HA - Adding a Secondary Appliance

Our recommended configuration is to use a clustered HA pair of load balancers to provide a highly available and resilient load balancing solution. We recommend that the Primary appliance is fully configured first, then the Secondary appliance can be added to create an HA pair. Once the HA pair is configured, load balanced services must be configured and modified on the Primary appliance. The Secondary appliance will be automatically kept in sync.

Note

For Enterprise Azure, the HA pair should be configured first. For more information, please refer to the Azure Quick Start/Configuration Guide available in the documentation library

The clustered HA pair uses Heartbeat to determine the state of the other appliance. Should the active device (normally the Primary) suffer a failure, the passive device (normally the Secondary) will take over.

13.1.1. Non-Replicated Settings

A number of settings are not replicated as part of the Primary/Secondary pairing process and therefore must be manually configured on the Secondary appliance. These are listed by WebUI menu option in the table below:

WebUI Main Menu Option	Sub Menu Option	Description
Local Configuration	Hostname & DNS	Hostname and DNS settings
Local Configuration	Network Interface Configuration	Interface IP addresses, bonding configuration and VLANs
Local Configuration	Routing	Default gateways and static routes
Local Configuration	System Date & time	Time and date related settings
Local Configuration	Physical – Advanced Configuration	Various appliance settings
Local Configuration	Portal Management	Portal management settings
Local Configuration	Security	Security settings
Local Configuration	SNMP Configuration	SNMP settings
Local Configuration	Graphing	Graphing settings
Local Configuration	License Key	Appliance licensing
Maintenance	Backup & Restore	Local XML backups
Maintenance	Software Updates	Appliance software updates
Maintenance	Firewall Script	Firewall (iptables) configuration
Maintenance	Firewall Lockdown Wizard	Appliance management lockdown settings

(!) Important

Make sure that where any of the above have been configured on the Primary appliance, they're also configured on the Secondary.


13.1.2. Configuring the HA Clustered Pair

8 Note

If you have already run the firewall lockdown wizard on either appliance, you'll need to ensure that it is temporarily disabled on both appliances whilst performing the pairing process.


- 1. Deploy a second appliance that will be the Secondary and configure initial network settings.
- 2. Using the WebUI on the Primary appliance, navigate to: *Cluster Configuration > High-Availability Configuration*.

Create a Clustered Pair

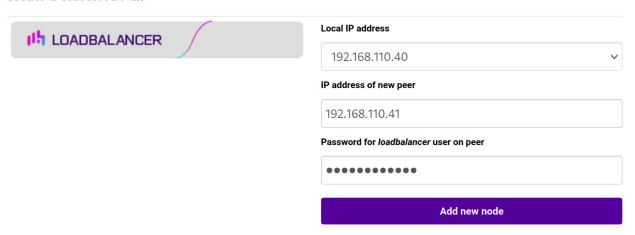

- 3. Specify the IP address and the *loadbalancer* user's password for the Secondary (peer) appliance as shown in the example above.
- 4. Click Add new node.
- 5. The pairing process now commences as shown below:

Create a Clustered Pair

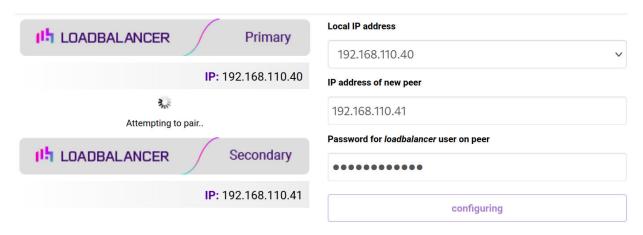

6. Once complete, the following will be displayed on the Primary appliance:

High Availability Configuration - primary



7. To finalize the configuration, restart heartbeat and any other services as prompted in the "Commit changes" message box at the top of the screen.



14. Document Revision History

Ve	ersion	Date	Change	Reason for Change	Changed By
1.	0.0	21 October 2025	Initial version		RJC



Visit us: www.loadbalancer.org

Phone us: +44 (0)330 380 1064

Phone us: +1 833 274 2566

Email us: info@loadbalancer.org

Follow us: @loadbalancer.org

About Loadbalancer.org

Loadbalancer.org's mission is to ensure that its clients' businesses are never interrupted. The load balancer experts ask the right questions to get to the heart of what matters, bringing a depth of understanding to each deployment. Experience enables Loadbalancer.org engineers to design less complex, unbreakable solutions and to provide exceptional personalized support.

